

ZA180A Zero air

Нулевой воздух

Номер части: 30-1660 (230v) - UK | 30-1670 (230v) - EU | 30-1680 (110v) - US




Описание

Генерируя нулевой воздух для детекторного газа в приложениях ГХ, генератор ZA180A удаляет все углеводороды из источника сжатого воздуха до того, как он попадет в ваше оборудование, обеспечивая самые чистые результаты от вашего детектора.



Generator Details

- **Типы газа:** Воздух
- **Встроенный компрессор:** да

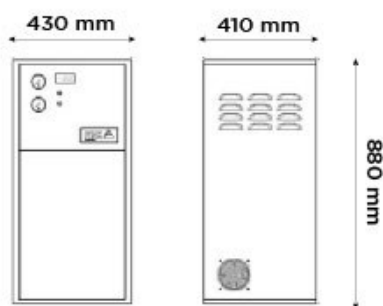

Ключевая особенность

- Производит нулевой воздух до содержания <math><0,1 \text{ ppm}</math>
- При необходимости работает круглосуточно и без выходных.
- Один генератор может питать несколько приборов.
- Меньшее количество движущихся частей означает минимальное техническое обслуживание и обеспечивает долгий срок службы генератора.
- Требуется минимальная настройка.
- 12 месяцев полной гарантии.

Технические характеристики

	ZA180A
Тип газа	Воздух
Макс. расход газа:	18L/min
Содержание углеводородов:	0.1ppm
Максимальная чистота:	0.1ppm THC
Время запуска:	40 Minutes
Потребляемая мощность:	770 Watts
Напряжение:	110 / 230 ± 10% VAC
Частота:	60 Hz
Текущий:	2/4 Amps
Максимальная рабочая температура:	30°C / 86°F
Аккредитации:	CE
Размер (ВхШхГ), мм	880 x 430 x 410 mm
Размер (ВхШхГ), дюймы	34,6 x 16,9 x 16,1 inches
Вес генератора	64kg / 140,8lbs

[PEAK Protected]TM

Peak Scientific gas generators define the benchmark in reliability, convenience and performance in laboratories around the world, and come backed with a 12 month warranty. Beyond this period however you can ensure that your investment continues to be **[Protected]** by our comprehensive generator care cover.

Our world-class aftercare support packages deliver a program of scheduled preventative maintenance whilst giving you the reassurance of instant access to worldwide technical support and priority on-site response in the untimely event of a breakdown.

- For ordering parts visit: www.peakscientific.com/ordering/
- For service plans visit: www.peakscientific.com/service/service-plans/

Peak Scientific UK
Tel: +44 (0)141 812 8100
Fax: +44 (0)141 812 8200

Peak Scientific North America
Tel: +1-800-767-6532
Fax: +1-978-608-9503

Peak Scientific China
Tel: +86 21 5079 1190
Fax: +86 21 5079 1191
For a full list of our worldwide office locations, please visit:
Web: www.peakscientific.com
Email: discover@peakscientific.com

Peak Scientific's Quality Management System conforms to: ISO:9001

